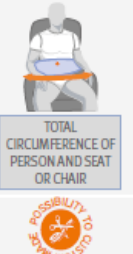


PRODUCT NAME:
ABDOMINAL BELT AND PERINEUM PIECE – ARNETEC® – Restraining Harnesses.
REFERENCE:
1005
SIZES:

SIZES	PERIMETER CM	
S/1	94-182	
L/2	104-192	

PRODUCT PICTURE:

MANUFACTURER:

ORLIMAN, S.L.U. CIF: B-96122510
 C/ AUSIAS MARCH Nº 3. POLÍGONO INDUSTRIAL LA POBLA – L'ELIANA.
 46185 LA POBLA DE VALLBONA – VALENCIA – SPAIN
 SPANISH MEDICAL DEVICE OPERATION LICENCE NUMBER: 20-PS

PRODUCT CLASSIFICATION (According to European Medical Device Regulation 2017/745):

CLASS: I – Annex VIII, Rule 1
 Manufacturing type: Mass-Produced.

PRODUCT CHARACTERISTICS:

Made from a triple layer of breathable semi-rigid material (polyamide, foam and polyester), it has regulation and securement straps with colour-coded buckle fasteners to facilitate positioning. The pelvic belt with perineum strap secures the patient to the wheelchair or other type of chair or seat, but at the same time allows full freedom of movement in the dorsal cervical zone.

Colour: Blue.

INDICATIONS:

- Suitable for patients with good anterolateral control whilst seated but with a tendency to slide out of the seat.

MATERIALS:

ALLERGEN: LATEX FREE
FERROMAGNETICS: NO
BIOCIDES: NO

Composition:

- Main fabric: 41% Polyester, 32% Polyurethane, 27% Polyamide
- Tape: 100% Polyester
- Fastening or catch elements:
 - Buckles: 100% POM
 - Straps: 100% Polyester

RECOMMENDATIONS FOR PRESERVATION AND CLEANING:

When not using the product, store it in the original packaging in a dry place at room temperature. Stick the Velcro to each other (if the orthotic device has them), frequently wash by hand with warm water (30° C max.) and mild soap. To dry the product, use a dry towel to absorb as much moisture as possible and let it dry at room temperature. Do not hang it up or iron the product and do not expose it to direct heat sources such as stoves, dryers, direct sun exposure, etc. When using or cleaning the product, do not use abrasive or corrosive substances, alcohol, ointments or liquid solvents. If not dried off properly, the detergent residue may irritate the skin and cause the product to deteriorate.

**ADDITIONAL INFORMATION:**

Guarantee: 6 months.

ACCESORIES**Ref.: 1100****PROLONGATION STRAP****Characteristics:**

Prolongation strap if needed for a wider seat, chair or bed (minimum 9 cm, maximum prolongation 85 cm). Except models magnets.

Colour: ■ Blue.**Packaging:** 1 unit.