

#040

PATELLA STABILIZER



CHARACTERISTICS AND FUNCTIONALITIES

Helps recover as fast and better as possible:

- Two removable patella c-shaped buttresses, that works in combination with the central straps, to optimize patella stabilization (medialization, centralization or lateralization)
- 3 elastic straps to adjust compression according to leg shape and recovery necessity
- Two removable lightweight and flexible lateral stays, to improve brace stability
- Free popliteal area - to avoid compression in this critical area and warranty normal blood circulation

User friendly:

- Ambidextrous

Comfortable:

- Breathable

Hygienic:

- 100% cotton in contact with the skin

INDICATIONS

- Patellar sub-luxation
- Patellar instability
- Tendonitis
- Post fractures
- Rehab passive support

#041

PATELLA AND MCL-LCL STABILIZER



CHARACTERISTICS AND FUNCTIONALITIES

Helps recover as fast and better as possible:

- Two removable patella c-shaped buttresses, that works in combination with the central straps, to optimize patella stabilization (medialization, centralization or lateralization)
- 3 elastic straps to adjust compression according to leg shape and recovery necessity
- Two removable polycentric shaped stays (one medial and one lateral) to give lateral stability to the knee and support MCL-LCL
- Free popliteal area - to avoid compression in this critical area and warranty normal blood circulation

User friendly:

- Ambidextrous

Comfortable:

- Breathable


Hygienic:

- 100% cotton in contact with the skin

INDICATIONS

- Patellar sub-luxation
- Patellar instability
- Tendonitis
- Post fractures
- Rehab passive support
- Light sprains
- Protective and contention action

Consisting of two C-shaped buttresses encircling the patella and stabilizing it.

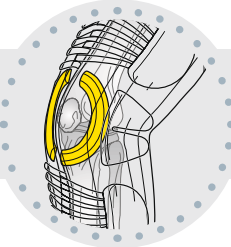

To exert an inward or outward pressure on the patella, just remove the lateral or medial buttress respectively.



AMBIDEXTROUS - HEIGHT: 31 CM

SIZE	S	M	L	XL	XXL	XXXL	XXXXL
(CM)*	32 - 36	36 - 40	40 - 45	44 - 49	48 - 54	54 - 60	60 - 66

* 10 cm Above The Knee

COLOURS

100% COTTON ON YOUR SKIN



100% Breathable



COLOURS

100% COTTON ON YOUR SKIN



100% Breathable

